HOMEDICS
THERAPEUTIC

Automatic Wrist Blood Pressure Monitor

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TABLE OF CONTENTS

INTRODUCTORY INFORMATION:

Important Product Notices and Safety Instructions3
About Blood Pressure6
What is Blood Pressure?	
Why Measure Your Blood Pressure?	
Blood Pressure Standard7
How This Blood Pressure Monitor Works8

Important Notes Regarding Your Blood Pressure Measurement27
Care, Maintenance & Cleaning29
Potential for Radio/Television Interference31
Troubleshooting33
Specifications34
Five Year Limited Warranty35

USAGE INFORMATION:

Name/Function of Each Part10
Display Explanations11
Installing Batteries13
Date & Time Set Procedure15
Applying the Cuff16
Correct Measurement Posture18
Measurement Procedure20
WHO Health Indicator24
Recalling Values from Memory25
Clearing Values from Memory26

IMPORTANT PRODUCT NOTICES AND SAFETY INSTRUCTIONS

When using your blood pressure monitor, basic precautions should always be followed. Please read and follow all instructions and warnings before using this product. Save these instructions for future reference.

- Please note that this is a home health-care product only and it is not intended to serve as a substitute for the advice of a physician or medical professional.
- Do not use this device for diagnosis or treatment of any health problem or disease. Measurement results are for reference only. Consult a healthcare professional for interpretation of pressure measurements. Contact your physician if you have or suspect any medical problem.

- This product is not suitable for people with arrhythmias. This device may have difficulty determining the proper blood pressure for individuals with diabetes, poor circulation of blood, kidney problems or for users who have suffered from a stroke.
- Consult your physician before measuring blood pressure at the wrist if you have any of the following conditions: Severe hypertension, severe diabetes, severe arteriosclerosis, severe kidney disease or vascular disease that may compromise circulation.
- The pulse display is not suitable for checking the frequency of heart pacemakers.
- Electromagnetic interference: The device contains sensitive electronic components. Avoid strong electrical or electromagnetic fields in the direct vicinity of the device (e.g. mobile telephones, microwave ovens). These may lead to temporary impairment of measurement accuracy.

- Use blood pressure monitor only for its intended use.
- Do not wrap the cuff around body parts other than your wrist.
- This unit is intended for use in measuring the blood pressure and pulse rate of adults. Do not use this device on infants or those who cannot express their own consent. This unit should not be operated by children.

Blood pressure measurements determined with this device are equivalent to those obtained by a trained observer using the cuff/stethoscope auscultatory method, within the limits prescribed by the American National Standard, *Manual, electronic, or automated sphygmomanometers*.

ABOUT BLOOD PRESSURE

What is Blood Pressure?

Blood pressure is the pressure exerted on the artery walls while blood flows through the arteries. The pressure measured when the heart contracts and sends blood out of the heart is systolic (highest) blood pressure. The pressure measured when the heart dilates with blood flowing back into the heart is called diastolic (lowest) blood pressure.

Why Measure Your Blood Pressure?

Among today's various health problems, those associated with high blood pressure are very common. High blood pressure dangerously correlates with cardiovascular diseases. Therefore, blood pressure monitoring is important for identifying those at risk.

BLOOD PRESSURE STANDARD

The **World Health Organization (WHO)** and the **National High Blood Pressure Education Program Coordinating Committee** have developed a blood pressure standard that identifies areas of low- and high-risk blood pressure. This standard, however, is a general guideline as blood pressure varies among different people, age groups, etc.

	Systolic (mmHg)	Diastolic (mmHg)
Stage 2 Hypertension	≥ 160	≥ 100
Stage 1 Hypertension	140~159	90~99
Prehypertension	120~139	80~89
Normal	< 120	< 80

To determine category of risk when systolic and diastolic readings fall into two areas, use the higher of the two numbers for classification.

This monitor comes equipped with a WHO Health Indicator to visually indicate the assumed risk level of each reading based upon these WHO standards. See page 24 for more information on this feature.

It is important that you consult with your physician regularly. Your physician will tell you your normal blood pressure range as well as the point at which you will be considered at risk. For reliable monitoring and reference of blood pressure, keeping long-term records is recommended.

HOW THIS BLOOD PRESSURE MONITOR WORKS

This monitor uses Smart Sense™ technology to detect your blood pressure. With one touch of a button, the cuff will automatically inflate to block the blood flow through your artery. Next, the deflation process starts. Smart Sense™ technology utilizes an advanced air control system to delicately control the deflation rate. The monitor efficiently analyzes and detects pulse signals while adjusting cuff pressure accordingly. Any muscle movement during this period of time will cause measurement error. When measurement is complete, the monitor will


display your systolic pressure, diastolic pressure, and pulse readings.

Based upon your measurement results, the monitor's WHO Health Indicator will alert you if your reading falls into the pre-hypertension, stage 1 hypertension, or stage 2 hypertension categories as classified by the World Health Organization. See page 24 for more information on the WHO Health Indicator.

NAME/FUNCTION OF EACH PART


Accessories for unit :

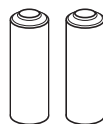2 "AAA" size, 1.5V alkaline batteries included in the carrying case with the unit.

DISPLAY EXPLANATIONS

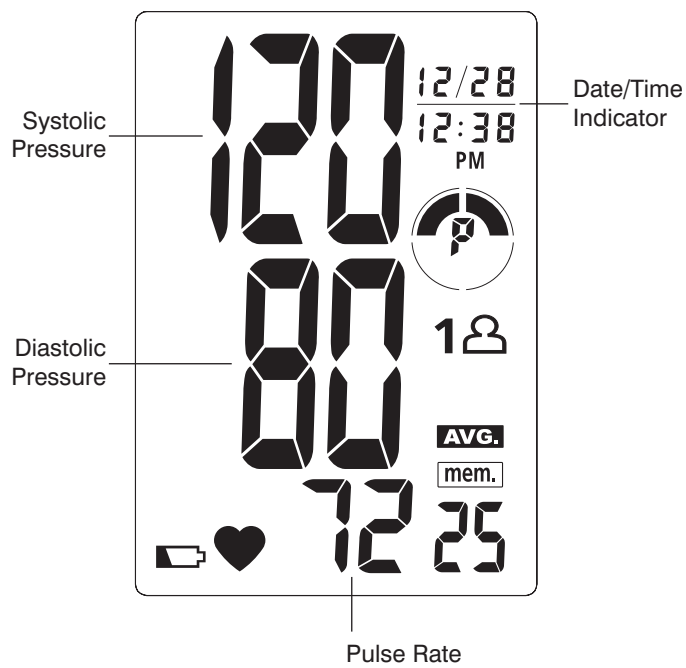Display Symbols:

1 (person icon)	User 1: Appears when the monitor is operated by User 1.
2 (person icon)	User 2: Appears when the monitor is operated by User 2.
(battery icon)	Weak Battery Symbol: Appears when batteries should be replaced.
(heart icon)	Pulse Symbol: Shows the pulse rate per minute.
mem.	Memory Symbol: Shows the memory reading number
AVG.	Memory Average: Displays average of last 3 readings
(WHO icon)	WHO Health Indicator: Compares blood pressure against guidelines set by the World Health Organization. See page 24 for more information.

If any of the following letters and numbers appear in the area that diastolic pressure should be displayed, an error has occurred with your reading. See Troubleshooting section of this manual for more information.





EE	Measurement Error: Measure again. Wrap the cuff correctly and keep wrist steady during measurement.
E1	Air Circuit Abnormality: Measure again.
E2	Pressure Exceeding 300 mmHg: Switch the unit off to clear, then measure again.
E3	Data Error: Remove and reload the batteries.

INSTALLING BATTERIES

1. The battery cover is located along the left side of the monitor. Place thumb and index finger on top and bottom ridges.
2. Remove battery cover by pressing and pulling cover away from the monitor.
3. Install lower battery first according to the polarity indications inside the compartment. This battery should be in the bottom position as shown in figure 3.


Battery Type: 2 Alkaline LR03 (AAA) size

4. Place in the other battery, according to the polarity indications inside the compartment and as shown in figure 4.

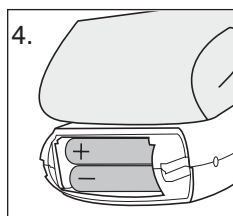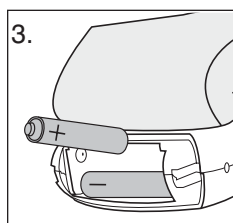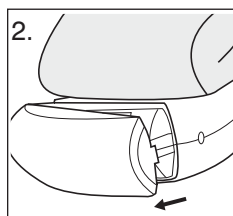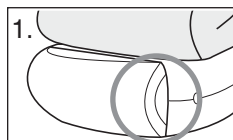13

5. Replace the cover by snapping it back into place.

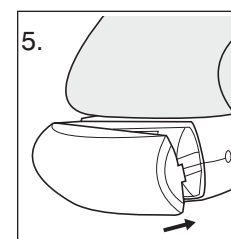Replace the batteries if:

1. The weak battery symbol appears in the display.
2. Nothing appears in the display when the power is switched on.


As the supplied batteries are for test only, they may be discharged earlier than batteries you buy in stores. Replace the batteries in pairs. Remove batteries when unit is not in use for extended periods of time.

Note: Batteries are hazardous waste. Do not dispose of them together with the household garbage.

14

DATE & TIME SET PROCEDURE


1. To adjust the date/time, press the **Set button** . The display will show a blinking number showing the month.
2. Change the month by pressing the **+** **button**. Each press will increase the number by one in a cycling manner. Press the **Set button**  again to confirm the entry, and the screen will show a blinking number representing the day of the month.
3. Change the date and time as described in Step 2 above, using the **+** **button** to change the numbers and the **Set button**  to confirm the entries.
4. A "0" will reappear as the blood pressure monitor is ready for measurement again.

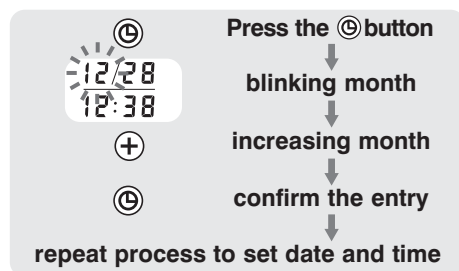15


APPLYING THE CUFF

1. Remove all watches, wrist jewelry, etc. prior to attaching the wrist monitor. Clothing sleeves should be rolled up and the cuff should be wrapped around bare skin for correct measurements.
2. Apply cuff to left wrist with palm facing up as shown in Fig. A.
3. Make sure the edge of the cuff is about 1/4"~1/2" (1 cm) from the palm as shown in Fig. B.

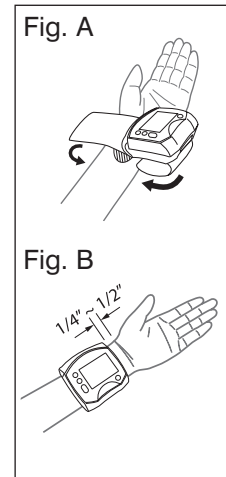
16

4. In order to insure accurate measurements, fasten the hook and loop strap securely around your wrist so there is no extra space between the cuff and the wrist as shown in Fig. C. If the cuff is not wrapped tight enough, the measurement values will not be accurate.

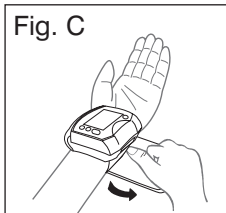5. If your physician has diagnosed you with poor circulation in your left arm, place the cuff around your right wrist as shown in Fig. D.


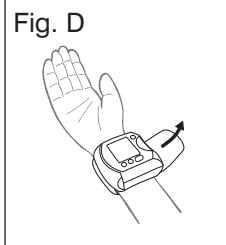
CORRECT MEASUREMENT POSTURE

1. Place your elbow on a table so that the cuff is at the same level as your heart as shown in Fig. A.

Note: Your heart is located slightly below your armpit.


Relax your entire body, especially between your elbow and fingers.

2. If the cuff is not at the same level as your heart or if you cannot keep your arm completely still throughout the reading, use a soft object such as a folded towel to support your arm as shown in Fig. B.

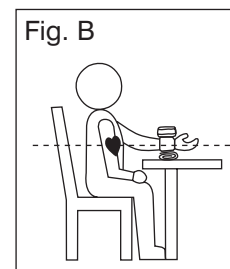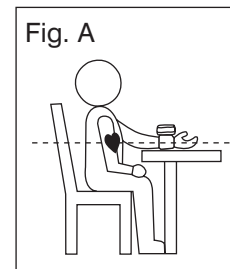3. Turn you palm upwards.
4. Sit upright in a chair, and take 5-6 deep breaths.

Avoid leaning back while the measurement is being taken as shown in Fig. C.

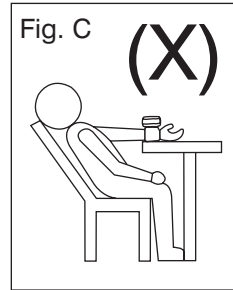
MEASUREMENT PROCEDURE

Important Notes:


- Position the monitor at the same level as your heart during measurement to ensure accurate readings.
- Blood pressure changes with every heart-beat and is in constant fluctuation throughout the day.
- Blood pressure measurement can be affected by the position of the user, his or her physiologic condition and other factors. For greatest accuracy, wait 1 hour after exercising, bathing, eating, drinking beverages with alcohol or caffeine, or smoking to measure blood pressure.
- Before measurement, it is suggested that you sit quietly for 15 minutes as measurements taken during a relaxed state will have greater accuracy. You should not be physically tired or exhausted while taking a measurement.
- During the measurement, do not talk or move your arm or hand muscles.
- See page 27 for additional notes regarding your blood pressure measurement

1. Press the **User Select Button** to choose User 1 or User 2.

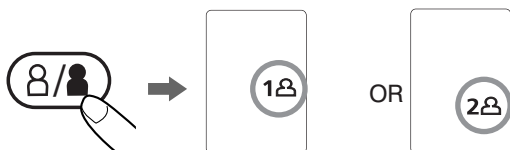
After the user number is selected, press the **START/STOP button** to confirm the chosen user.

2. With the cuff wrapped around your wrist, press the **START/STOP button**. Do not inflate the cuff unless it is wrapped around your wrist. All digits will light up, checking the display functions. The checking procedure will be completed after about 2 seconds.

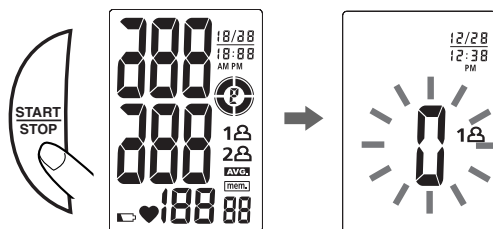
21

3. After all symbols disappear, the display will show a blinking "0". The monitor is **"Ready to Measure"** and will **automatically** inflate the cuff to begin measurement. The cuff will then begin deflating as measurement continues.

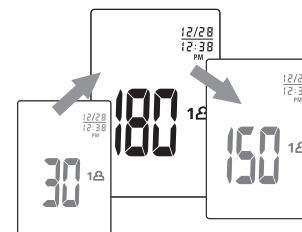
Note: The monitor will re-inflate automatically if the system detects that your body requires more pressure for measurement.

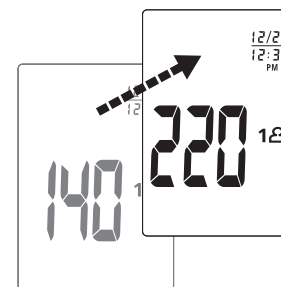22

4. When the measurement is completed, the cuff will deflate entirely and systolic pressure, diastolic pressure and pulse will be shown simultaneously on the LCD screen. **The measurement is then automatically stored into memory.**

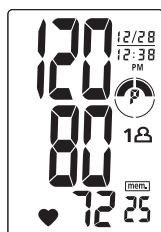Note:

1. This monitor automatically turns off approximately **1 minute** after last operation. You may also press the **START/STOP button** to turn the unit off.
2. To interrupt the measurement, you may press the **START/STOP button** (recommended) or the Memory, User Select, or Date/Time Set (Ⓢ) buttons. The cuff will deflate immediately after a button is pressed.

WHO HEALTH INDICATOR

The World Health Organization (WHO) and the **National High Blood Pressure Education Program Coordinating Committee** have developed a blood pressure standard, classifying blood pressure ranges into 4 stages. This unit is equipped with an innovative blood pressure classification indicator, which visually indicates the assumed risk level (normal, prehypertension, stage 1 hypertension, or stage 2 hypertension) of each reading. Symbols are described below and classifications follow the WHO guidelines as described earlier in this manual (page 7).

	Stage 2 Hypertension		Pre-hypertension
	Stage 1 Hypertension		Normal (No Symbol)

RECALLING VALUES FROM MEMORY

This monitor can be used by two individuals.
Each user can store up to 60 measurements.

1. Press the **User-Select button**.

2. Choose User 1 or User 2.


3. Press the **“M” button** to access the memory.

4. **The monitor will first display the calculated average applied to the last 3 memories.**

5. Every new press of the **“M” button** will recall a previous reading. The latest reading will be recalled first.

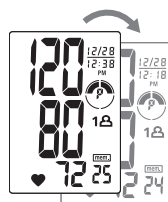
CLEARING VALUES FROM MEMORY

1. Press the **User-Select button** to select User 1 or User 2.

2. Continue to press the **“M” button** for approximately 5 seconds, then the data for the selected user will be erased automatically.

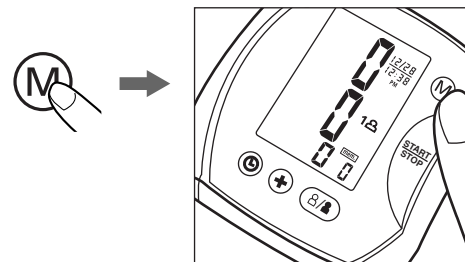Note: The data stored in memory for both users may be erased if any of the batteries are removed. (Users should first record the data in a Blood Pressure log before replacing or removing batteries.)

IMPORTANT NOTES REGARDING YOUR BLOOD PRESSURE MEASUREMENT

- Take your reading in a comfortable environment as measurements can be affected by hot or cold temperatures. Take your blood pressure at normal body temperature.
- Do not move or talk during measurement as this can elevate readings.
- Do not move or cross legs during measurement. Keep feet flat on floor.
- Do not touch cuff or monitor during measurement procedure.
- It is suggested that you take your measurements at the same time each day and use the same arm for consistency.
- Users should wait a minimum of 5 minutes before taking additional measurements. More time may be necessary depending upon your physiology.
- **The measurement results that users receive are for reference only. If users have any blood pressure concerns, please consult a physician.**

- Once inflation reaches 300 mmHg, the unit will deflate automatically for safety reasons.
- **This product is not suitable for people with arrhythmias.**
- **This device may have difficulty determining the proper blood pressure for users with diabetes, poor circulation of blood, kidney problems, or for users who have suffered a stroke.**

CARE, MAINTENANCE & CLEANING

1. Clean the blood pressure monitor body and cuff carefully with a slightly damp, soft cloth. Do not press. Do not wash cuff or use chemical cleaner on it. Never use thinner, alcohol or petrol (gasoline) as cleaner.
2. Leaky batteries can damage the unit. Remove the batteries when the unit will not be used for a long time.
3. If the unit is stored near freezing, allow it to acclimate to room temperature before use.
4. This blood pressure monitor is not field serviceable. You should not use any tool to open the device nor should you attempt to adjust anything inside the device. If you have any problems with this device, please contact HoMedics Consumer Relations (contact information on warranty page).
5. Do not immerse the unit in water as this will result in damage to the unit.

6. Do not subject the monitor and cuff to extreme temperatures, humidity, moisture, and direct sunlight. Protect from dust.
7. Do not fold the cuff tightly.
8. Do not disassemble the monitor or cuff. If in need of repair, refer to the warranty section of this manual.
9. Do not subject the monitor to extreme shocks (do not drop on floor).
10. Do not inflate the cuff unless wrapped around wrist.
11. Do not wrap the cuff around body parts other than your wrist.
12. Do not drop or insert any object into any opening.
13. This monitor may not meet its performance specifications if stored or used outside of these temperature and humidity ranges:

• Storage Environment	• Operation Environment
Temperature: 14°F~140°F (-10°C~60°C)	Temperature: 50°F~104°F (10°C ~40°C)
Humidity: 10~90% RH max.	Humidity: 40~85% RH max.

POTENTIAL FOR RADIO/TELEVISION INTERFERENCE

For USA only

This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of FCC rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. The product generates, uses, and can radiate radio frequency energy and, if not used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If the product does cause harmful interference to radio or television reception, which can be determined by turning the product on and off, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the product and the receiver.
- Consult the dealer or an experienced radio/TV technician for help.

POTENTIAL FOR RADIO/TELEVISION INTERFERENCE

For Canada only

This digital apparatus does not exceed the Class B limits for radio noise emissions from digital apparatus as set out in the interference-causing equipment standard entitled "Digital Apparatus", ICES-003 of the Canadian Department of Communications.

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TROUBLESHOOTING

If any abnormality arises during use, please check the following points:

SYMPTOMS	POSSIBLE CAUSES	CORRECTION
Unit does not turn on when the START/STOP button is pushed.	Batteries have run down.	Replace them with two new AAA alkaline batteries.
	Battery polarities have been positioned incorrectly.	Re-insert the batteries in the correct positions.
EE measurement error symbol shown on display or the blood pressure value is displayed excessively low (or high).	Cuff has been placed incorrectly.	Wrap the cuff properly so that it is positioned correctly.
	Did you talk or move during measurement?	Measure again. Keep wrist steady during measurement.
	Shaking of the wrist with the cuff on.	
E1 error symbol shown on display	Air circuit abnormality.	Measure again.
E2 error symbol shown on display	Inflation pressure exceeding 300 mmHg.	Switch the unit off, then measure again.
E3 error symbol shown on display	Data error.	Remove and reload the batteries.

Note: If the unit still does not work, contact HoMedics Consumer Relations. Under no circumstance should you disassemble or attempt to repair the unit by yourself. Contact information for HoMedics Consumer Relations Department can be found on the warranty page (See page 35).

SPECIFICATIONS

Power Source:	Two 3V DC (AAA) Alkaline batteries
Measurement Method:	Oscillometric
Measurement Range:	Pressure: 40~250 millimeters Mercury (mmHg) Pulse: 40~199 beats/minute
Accuracy:	Pressure: ± 3 mmHg Pulse: $\pm 5\%$ of reading
Pressure Sensor:	Semi-conductor
Inflation:	Pump driven
Deflation:	Automatic Pressure release valve
Memory Capacity:	60 memories for each User (120 total)
Auto-shut-off:	1 minute after last button operation
Operation Environment:	Temperature: 50°F~104°F (10°C~40°C) Humidity: 40%~85% RH max.
Storage Environment:	Temperature: 14°F~140°F (-10°C~60°C) Humidity: 10%~90% RH max.
Weight:	0.30 lbs (without batteries)
Wrist Circumference:	5.3"~8.7" (13.5~22cm)
Dimensions:	2.6"(L) x 1.3"(W) x 3.1"(H) 67mm(L) x 32mm(W) x 78mm(H)
	Type BF Device and cuff are designed to provide special protection against electrical shocks.

Note: These specifications are subject to change without notice.

HoMEDICS
THERAPY

For service or repair, do not return this unit to the retailer. Contact HoMedics Consumer Relations at:

Mail:
HoMedics, Inc.
Service Center
Dept. 168
3000 Pontiac Trail
Commerce
Township, MI
48390

Email:
cservice@
homedics.com

Phone:
1-800-466-3342
Business Hours:
8:30am-5pm ET
Monday-Friday

FIVE YEAR LIMITED WARRANTY (Valid in USA only)

HoMedics, Inc., guarantees this product free from defects in material and workmanship for a period of five years from the date of original purchase, except as noted below.

This HoMedics product warranty does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorized accessory; alteration to the product; or any other conditions whatsoever that are beyond the control of HoMedics. This warranty is effective only if the product is purchased and operated in the USA. A product that requires modification or adaptation to enable it to operate in any country other than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under warranty. HoMedics shall not be responsible for any type of incidental, consequential or special damages.

All implied warranties, including but not limited to those implied warranties of fitness and merchantability, are limited in the total duration of five years from the original purchase date.

5 year
limited warranty

HoMEDICS
THERAPY

To obtain warranty service on your HoMedics product, either hand deliver or mail the unit and your dated sales receipt (as proof of purchase), postpaid, along with check or money order in the amount of \$5.00 payable to HoMedics, Inc. to cover handling.

Upon receipt, HoMedics will repair or replace, as appropriate, your product and return it to you, postpaid. If it is appropriate to replace your product, HoMedics will replace the product with the same product or a comparable product at HoMedics' option. Warranty is solely through HoMedics Service Center. Service of this product by anyone other than HoMedics Service Center voids warranty.

This warranty provides you with specific legal rights. You may have additional rights which may vary from state to state. Because of individual state regulations, some of the above limitations and exclusions may not apply to you.

For more information regarding our product line in the USA, please visit: www.homedics.com

©2005 HoMedics, Inc. and its affiliated companies, all rights reserved. HoMedics is a registered trademark of HoMedics, Inc. and its affiliated companies. Smart Sense™ and TheraP™ are trademarks of HoMedics Inc. and its affiliated companies. All rights reserved.

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